

HWB041-10

PATIENT TRANSPORTATION TROLLEY



Healthward

Healthier your life

HWB041-10

FEATURES AND HIGHLIGHTS

Technical parameters

External size (LxW)	1880x620 mm
Table height adjustment	560-890 mm
Back plate adjustment	0-75°
Knee-rest adjustment	0-40°
Revrese trendelenburg	0-18°
Trendelenburg	0-18°

Technical configuration

1 set
1 set
1 set
4 pcs
1 pc
2 pcs
1 set
4 pcs
1 pc
1 pc





















Guardrail

Side-rail can be Fixed at Horizontal Level Fixable at horizontal level, where patient's arm can be placed on for the administration of IV or other treatment.

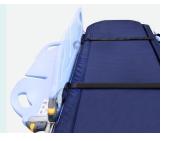

One direction moving

Proper procedure needs to be exercised to safely transfer the patient from the stretcher to the bed. An innovative side-rail design allows medical personnel to carry out this action efficiently and with greater ease. Side-rail can be laid out on a bed positioned next to the stretcher to act as a transfer board, creating zero gap for easy and quick transfer.

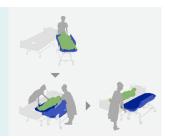

Degree indicater

Scale on guardrail, shows the backrest degree.



T-shaped gap

T-shaped design allows tubes to pass through side-rail for improved fluid drainage.



Hydraulic system

American brand hydraulic pump, providing strong support for the functions of height adjustment, trendelenburg/reversed trendelenburg.



Bumper

Equipped with bumper at four



Foldable handle

Equipped with foldable handle.



Center Wheel

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.



Caster

8" centrolled castors.



I.V pole

Full SS IV pole, bearing 15kg.



Mattress

Wator proof mattress.



Healthward

As the circumstance of printing, the photos might have minor difference on color, tactile from real products. As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.















export@health-ward.com

www.health-ward.com

