




HWB041-6

Patient
Transportation
Trolley



Healthward

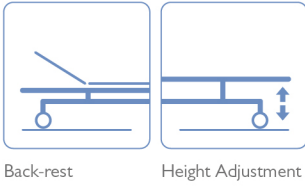

Back-rest

Height Adjustment

MAIN TECHNICAL PARAMETERS

Size	1930*660*(510-850)mm
Back-rest	0-70°(±10°)

TECHNICAL CONFIGURATION

Crank System	1 set
Back Rest Air Spring	1 pc
Plastic Side Rail	1 set
6" Castors , Central Locking	4 pcs
Directional Wheel	1 pc
Directional Pedal	2 pcs
IV Pole Previsions	2 pcs
IV Pole	1 set
Oxygen Cylinder Holder	1 pc
Central Locking Break Pedal	4 pcs
Transfer Mattress	1 pc
Safe Lock for Side Rail	2 pcs

QUALITY CONTROL

GUARDRAIL REVERSIBLE

Side-rail can be Fixed Horizontal Level. Fixable at horizontal level, where patient's arm can be placed on for the administration of IV or other treatment.



T SHAPE DESIGN

T-shaped design allows tubes to pass through side-rail for improved fluid drainage.



DEGREE INDICATE

Scale on guardrail, shows the backrest degree.



ONE DIRECTION MOVING

The side-rail can be laid out one a bed next to the stretcher to act as a transfer board, allows medical personnel transfer patient from stretcher to bed more efficiently and easier.




CENTER WHEEL

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.

RICH PARTS

LEADING BRAND OF
CHINA MEDICAL FURNITURE

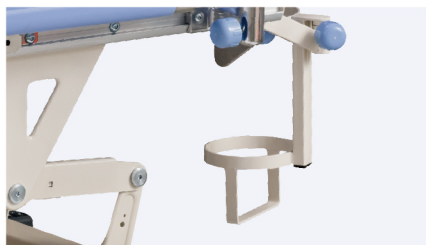

Adjustable vertical oxygen bottle, load capacity: 15KG



Stainless steel handle with plastic handle, more durable and beautiful



6 inch controlled castors



Full SS IV pole, bearing 15kg



Water proof mattress with 2 belts



As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.

Healthward



 17800 CASTLETON ST STE 180.CITY OF INDUSTRY CA 91748

 +1 (603) 416-3063 ·  24/7

 export@health-ward.com

 www.health-ward.com



General catalogue