

HWR-R02 MOBILE MULTI-FUNCTION INFORMATION SEARCH CART





HWR-R02

FEATURES AND HIGHLIGHTS

Technical Parameters

 Dimensions L×W×H 	620×540×1460mm
 Work table L×W 	460×480mm
Weight of the whole machine	45KG (without computer)
 Sub-table size L×W 	260×260mm
Trolley lift stroke	200mm
Display lifting stroke	90mm
Caster design	4 inch medical grade silent wheel
Cabinet parameters	two layers
Drawer dimensions	480×340×295mm
Layer drawer inner size	390×275×55mm
• The inner size of the drawer on the s	second floor 390×275×140mm

Optional Accessories

Hand sanitizer box
☐ Mouse hanger
Sundries box
Trash can



Material

The table top is made of ABS medical grade plastic, the two-layer allaluminum box is made of high-voltage electrostatic spraying process, and imported environmentally friendly antibacterial plastic powder is used to ensure that the product meets the medical and health requirements, and the surface is matte and non-dazzling. And equipped with an integrated die-casting aluminum handle and base,

The column material is made of medical grade aviation aluminum alloy.

Upper Column

The detachable upper column can carry a monitor or an all-in-one computer up to 10kg, and the universal support function can complete 180° left and right rotation and 45° up and down inclination. The up and down travel of the column is as high as 15cm, with a constant force lifting adjustment for easy operation, and various equipment threading can be hidden in the column to make the appearance more refreshing.

Table Top

The ABS medical-grade plastic countertop is easy to clean and anticorrosion, without any spraying on the surface, which will not affect the medicine. The grooves are set around to prevent the dripping of medical liquid, and there is a countertop enclosure to prevent items from falling.

Sub Table Top

The pull-out sub-table provides the area expansion of the desktop, which can be temporarily convenient for work.







Keyboard Drawer

Hidden keyboard drawer that can be stored inside the countertop when not in use.

Electromechanical Layer

The distance between the electromechanical layer and the medicine box is greater than 5cm, which facilitates the heat dissipation of the electromechanical equipment and protects the safety of medicines.

Drawer

Two-layer large space self-priming drawer, the inner layer of the drawer is equipped with a plastic insert board, which can be combined arbitrarily to store medicines of different specifications.


Handle

The ring-shaped anti-collision handle, surrounded by solid metal handles with rounded corners, is located at the outermost periphery of the whole vehicle, which can not only injure other personnel in the event of a collision, but also protect the safety of the stroller.







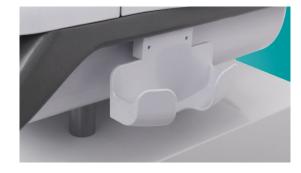

Lower Column

The column is equipped with a professionally adjusted air rod, which is easy to rise and fall, and can be quickly adjusted to a suitable height, adapting to different sitting and standing postures, with a stroke of up to 20cm. And optional electric lifting system.



Chassis

All-metal stretched X-shaped chassis, one-time cast aluminum molding, to ensure the center of gravity and level of the vehicle.



Casters

4-inch dual-wheel medical grade silent wheels, all with brakes, after tens of thousands of kilometers of long-term durability tests, flexible and durable.

Monitor (optional)

Split computer: 21.5 inches or 23.8 inches All-in-one computer: 23.8 inches and above





Healthward `

As the circumstance of printing, the photos might have minor difference on color, tactile from real products. As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.



- 17800 CASTLETON ST STE 180.CITY OF INDUSTRY CA 91748
- 🕕 +1 (603) 416-3063 · 🕓 24/7
- export@health-ward.com
- () www.health-ward.com

