


HWV3v MANUAL BED

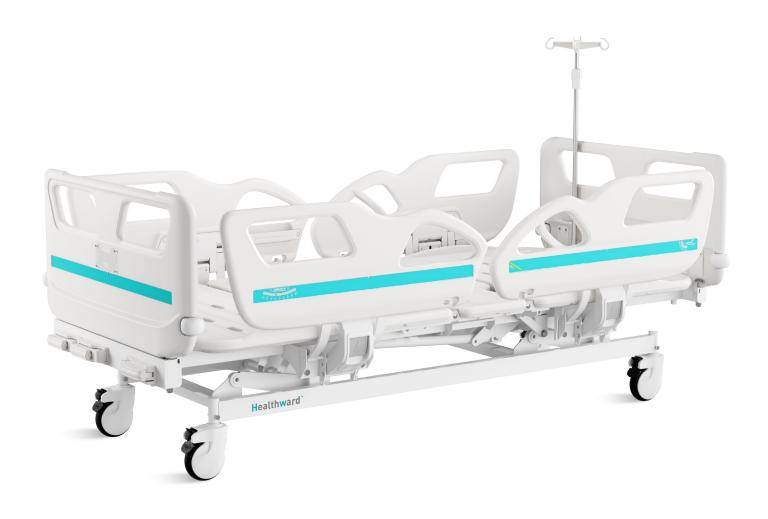

Healthward`

Healthier your life

HWV3v

FEATURES AND HIGHLIGHTS











Technical parameters

•	External size (LxW):	2130x1025 mm
	Mattress platform:	1925x900 mm
	Safe working load:	200 kg
	 Height adjustment: 	430-740 mm
•	■ Back-rest adjustment:	0-73° (±5°)
	Knee-rest adjustment:	0-38° (±5°)
	Angle indicator:	0-90°
•	• Auto-regression:	100 mm
•	Position line:	30°

Technical configuration

☑ 5" Covered Castors 4 pc	S
✓ PP side rail extra height 1 se	et
✓ PP bed ends with safe lock with Φ65 bumper wheels 1 se	et
✓ PP bed platform with integrated mattress retainer 1 se	et
✓ IV pole 4 hooks 1 see	et
☑ IV pole prevision 4 po	S
☑ Drainage hook 2 po	S

Optional

☐ Utility shelf

☐ Trapeze bracket/pole

Optional color (MOQ100pcs-color changing)



QUALITY CONTROL

Healthward, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. Healthward adopt robot welding, automatic epoxy coating. efficient assembly line, dedicated operator, rigorous quality control, only to guarantee Healthward's commitment on quality



Design

The unique germ-free design meets hospitals' strict infection control requirment. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization



Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Auto-Regression

Auto-Regression 10cm.
When repositioning the backrest and thigh rest, the V bed makes the back plate retreat 10cm, so to minimises the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.

Accessories



Platform

Integrated PP bed platform, adopting advanced manufacturing technology. 4-5 section platform, enables independent installation and easy maintenance.

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



Castors

5" covered casters. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.



















Back-rest adjustment



Knee-rest adjustment

38°



Height adjustment

 $430\text{-}740_{\mathsf{mm}}$













🗍 +1 (603) 416-3063 · 🕓 24/7

export@health-ward.com

www.health-ward.com

